


Rx, OTC, Nutraceutical

LAB SCALE SOFT GELATIN CAPSULE MACHINE

SGL107

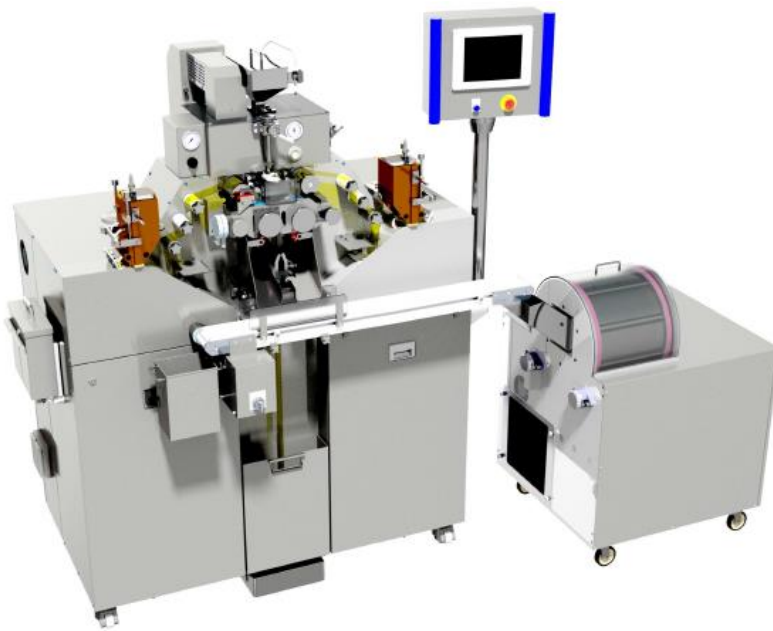

Rx, OTC, NutraceuticalSGL107

LAB SCALE SOFT GELATIN CAPSULE MACHINE



Features

- Laboratory scale soft capsule machine
- Motion Control System: the most important mechanical gearing replaced with electronic gearing in synchronized loops, using servo motors
- Die Roll dimensions: approximately 4" x 3"
- Operation speed - up to 3 RPM.
- Control of all production and process parameters using Programmable Logic Controllers (PLC) with touch screen
- Electrical panel and PLC incorporated within machine
- Sealed construction to prevent leakage and subsequent contamination
- Aluminum frame construction prevents any possible rust and allows easy cleaning and maintenance. The outside covers of the machine are completely made from stainless steel, designed according to the GMP/FDA requirements



Mechanism

Medicine Pump System

- The medicine pump is independently driven by a servo motor and electronically synchronized with the Die Rolls
- During a capsule size change the operator inputs the number of die roll cavities using the touch screen and the medicine pump timing is done electronically (no gear change is required)
- Precise injection timing adjustments are done electronically using the touch screen; The timing can also be adjusted during operation, without stopping the machine. Since this timing can be adjusted with such high accuracy, this improves on seal integrity of the capsule, which will increase capsule quality
- The slide lubrication system for the pump is integrated in the drive mechanism
- Single plunger standard medicine pump
- Single plunger Maxi and Mini medicine pumps available

Spreader Box

- Adaptable gelatin ribbon width that allows a minimum loss of ribbon net;
- Bronze construction
- Gelatin level control is done electronically for gravity feeding from individual attached tanks; ground level gelatin feeding also available, using a pressurized gelatin tank (optional)

Casting Drums

- “In-line” construction allows a straight gelatin ribbon transfer to the oil rolls and die rolls without any twisting
- 10.24" (260mm) diameter air cooled drums
- Water drum cooling system available

Ribbon Lubrication & Drive System

- Porous material for lubrication roller
- Electronic metering pump and a precise distributor valve allow accurate ribbon lubrication
- Two oil formulations and 4 metering pumps available
- Easy and accurate oil quantity adjustment
- Independent electrical drives allow fine adjustments of the ribbon tension on each side
- Easy access from the operator side to feed the ribbon on the oil rollers

Rx, OTC, Nutraceutical^{SG}L107

LAB SCALE SOFT GELATIN CAPSULE MACHINE



Mechanism cont'd

Die Roll Pressure System

- Manually operated pneumatic valves for High and low pressure

The Divider (Chute Assembly)

- Variable speed independent drive
- Sealed construction

Mangle Rolls

- Electronic speed adjustment allows fine tuning of the ribbon tension and helps in the cutting process

Wedge Lifting System

- Manually operated with safe middle position of the injection wedge;

Shut off valve

- 2 positions manual valve

Machine Lubrication

- Closed lubrication system with dedicated pump for the Die roll housing
- Separate lubrication pump for the medicine pump slide

Soft Capsule Transfer Conveyor

- Variable speed drive motor
- Bi-directional

Set of aluminum die rolls (4" x 3") and change parts


- Specially coated injection wedge
- Storage case
- Die size – per customer requirements

Rx, OTC, Nutraceutical^{SGL107}

LAB SCALE SOFT GELATIN CAPSULE MACHINE



Specification

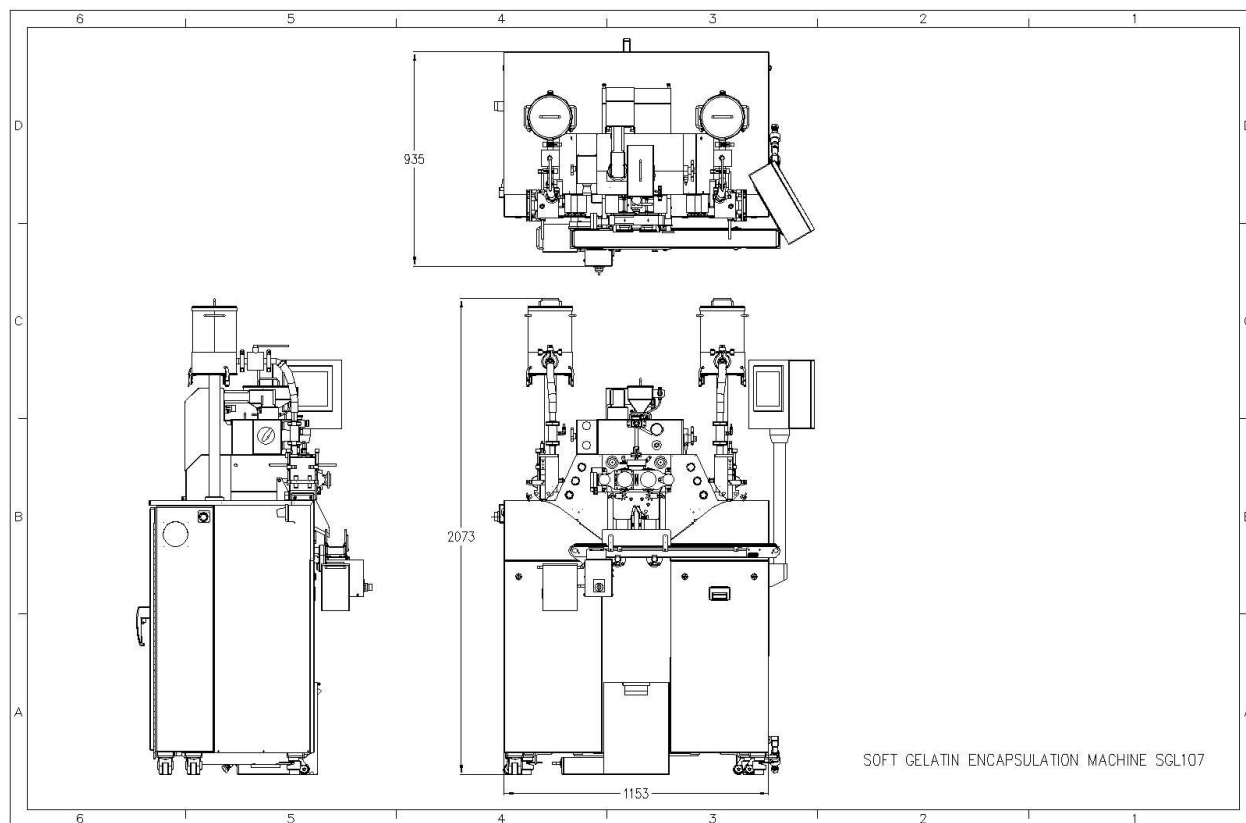

SGL107	
Dimensions (mm)	1153 W X 845 D X 2073 H
Weight (kg)	850
Machine Electrical Power	12 KVA
Compressed Air	100 PSI
Chilled Water	3.8 GPM

TECHNOPHAR

Technophar Global Headquarters
3293 St. Etienne Blvd.
Windsor, Ontario, Canada N8W 5B1
Phone: +1 519-946-0800
www.technophar.com